

SILICONE BIRTH-VAC™ CUP

Silicone Obstetric Vacuum Extractor A
safe and handy birth assistant



GYNECOLOGIC &
OBSTETRIC
SERIES

SILICONE BIRTH-VAC™ CUP

Silicone Birth-Vac™ Cup is designed for birth assisting. Made of silicone material, the cup comes with a soft rim which is made to minimize the risk of trauma to fetal head. A reference-line imprinted vertically on the cup for visual positioning. A rigid shaft with three loops provides firm grip. A chromed plated metal handle offers easy-hold, and the valve of the handle provides easy accessibility for vacuum adjustment.


Silicone Birth-Vac™ Cup shortens and assists second stage of labor in occipital delivery, which avoids the danger of forceps delivery. It also minimizes the risk of trauma and cephalohematoma of fetal cranium. Designed for reusable, Birth-Vac™ Cup is made of Medical grade, which features prolonged product life with proper use, cleaning and reprocessing procedures. Three cup sizes of 50, 60 and 70 mm are available for selection



SPECIFICATIONS

Birth-Vac™ Cup:

REF. No.	Size.	Application
1300-0050	50 mm	- Silicone
1300-0060	60 mm	- Soft rim
1300-0070	70 mm	- Rigid shaft

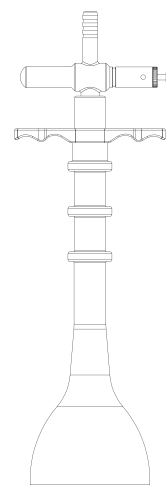


1300-00XX

GYNECOLOGIC & OBSTETRIC SERIES

SILICONE BIRTH-VAC™ CUP

Replacement Valve:

REF. No.	Size.	Application
1300-9001	Universal	Chrome plated metal

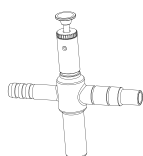

1300-9001

Connecting Tubing:

REF. No.	Size.	Application
1300-9200	6.0 x 11.0 mm	200 cm in length



1300-9200



6 F., No. 29, Sec. 2, Jhongjheng E. Rd.,
Danshuei Dist, New Taipei City. 251 Taiwan
Tel.:(886) 2-2624-2233 Fax.:(886) 2-2624-2266
www.fortunemed.com Mail:info@fortunemed.com